


Learn how Newman Medical started and developed the Natus Doppler line.



DigiDop

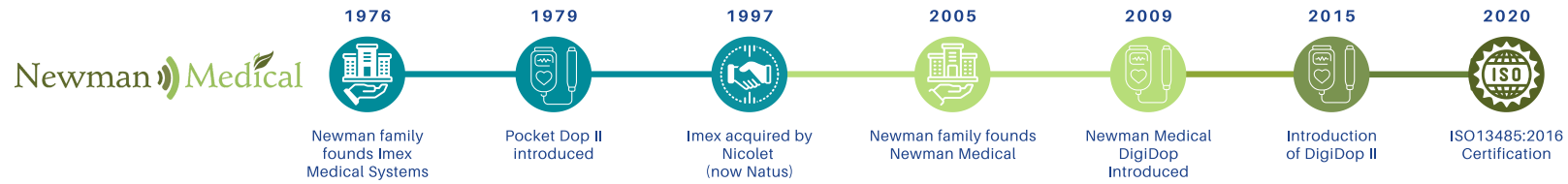

Includes labor, parts, and accidental breakage



American-Made Quality



Includes cords and probes



Waterproof Option



Rechargeable Option

H A N D H E L D D O P L E R S  
N O N - D I S P L A Y

natus Model



Elite 100

Elite 100/100R



PD-II

Pocket-Dop II



Elite 200

Elite 200/200R

Elite 200R



IMEXDOP CT+

CT+

REPLACEMENT  
Newman Medical  
DigiDop Classic



DD-300



DD-301



DD-301



DD-700



DD-701

REPLACEMENT  
Newman Medical Model



DD-901

REPLACEMENT  
Newman Medical  
DigiDop II



DD-330



DD-330R



DD-330R



DD-770



DD-770R

P R O B E O R D E R I N G I N F O R M A T I O N O N B A C K S I D E

## Vascular Probes



**Model D5**  
5MHz Vascular Probe



**Model D8**  
8MHz Vascular Probe



**Model DPPG**  
Audio TBI PPG Probe

## Fetal Probes



**Model D2**  
2MHz Fetal Probe



**Model D2W**  
2MHz Waterproof Fetal Probe



**Model D3**  
3MHz Fetal Probe



**Model D3W**  
3MHz Waterproof Fetal Probe

## Professional and Sensitive Dopplers

### Quality

- Better sound - Digital Sound Clarity™ produces up to 6x signal of other Dopplers
- Better warranty - 5 year honest warranty all parts & labor. Includes rechargeable batteries on rechargeable units.
- Better support & service - knowledgeable staff and guaranteed service

### Sustainability

- Smart Recharge System™ (SRS)
  - Significant cost savings
  - New rechargeable batteries hold charge for years
  - May be recharged more than 2000 times

### Innovation

- Continual product development and research

### DigiDop Classic



- Standard size
- Longer battery life than similar Natus products

### DigiDop II



- Smaller size
- Ergonomic - fits comfortably in your hand
- Similar battery life to Natus

### DigiDop 901



- Ideal for hospital use
- Great room-filling sound from two speakers on top
- Two probe cradles
- Standard drop-on charging pyramid
- Optional wall mount recharger

## How to build your Doppler

DOPPLER CODE

PROBE CODE

### Doppler

|        |                                  |
|--------|----------------------------------|
| DD-300 | Basic                            |
| DD-301 | Basic + Recharge                 |
| DD-700 | HR Display                       |
| DD-701 | HR Display + Recharge            |
| DD-901 | Tabletop + HR Display + Recharge |

### Probe

|     |                             |
|-----|-----------------------------|
| D2  | 2MHz fetal probe            |
| D2W | 2MHz waterproof fetal probe |
| D3  | 3MHz fetal probe            |
| D3W | 3MHz waterproof fetal probe |

DigiDop "Classic"

DigiDop II

|         |                       |
|---------|-----------------------|
| DD-330  | Basic                 |
| DD-330R | Basic + Recharge      |
| DD-770  | HR Display            |
| DD-770R | HR Display + Recharge |

|      |                                     |
|------|-------------------------------------|
| D5   | 5MHz vascular probe                 |
| D8   | 8MHz vascular probe                 |
| DPPG | Audio PPG probe - includes toe cuff |

## Interchangeable probes can be used on any DigiDop

### Vascular



- Great sensitivity
- Ergonomic shape easy to hold, place in sterile glove finger for use in sterile fields
- D8 - 8 MHz** - best for shallow vessels like radial, ankle, foot, digits, penile
- D5 - 5 MHz** - wider beam, deeper penetration. Good for brachial, carotid, femoral, deep veins

### Fetal



- Great sensitivity
- Smooth shape, comfortable for patient, easy to hold and clean
- D3 - 3 MHz** - wide angle, ideal for early fetal hearts
- D3W** - a unique 3 MHz waterproof probe
- D2 - 2MHz** - deeper penetration
- D2W** - 2 MHz - waterproof

### PPG



- DPPG**
- Unique audio PPG probe for use on digits. Great for toe-brachial studies

## NATUS-NEWMAN MEDICAL PART NUMBER CROSS-REFERENCE



**Elite**  
**100/100R**

| Natus | DigiDop Classic | DigiDop II |
|-------|-----------------|------------|
| EN20  | DD-300-D2       | DD-330-D2  |
| EN20R | DD-301-D2       | DD-330R-D2 |
| EN30  | DD-300-D3       | DD-330-D3  |
| EN30R | DD-301-D3       | DD-330R-D3 |
| EN50  | DD-300-D5       | DD-330-D5  |
| EN50R | DD-301-D5       | DD-330R-D5 |
| EN58  | DD-300-VASC     | --         |
| EN80  | DD-300-D8       | DD-330-D8  |
| EN80R | DD-301-D8       | DD-330R-D8 |
| EN8A  | DD-PAD          | --         |



**Elite**  
**200/200R**

| Natus | DigiDop Classic | DigiDop II  |
|-------|-----------------|-------------|
| ED20  | DD-700-D2       | DD-770-D2   |
| ED20R | DD-701-D2       | DD-770R-D2  |
| ED23  | DD-700-OB       | --          |
| ED2W  | DD-700-D2W      | DD-770-D2W  |
| ED2WR | DD-701-D2W      | DD-770R-D2W |
| ED30  | DD-700-D3       | DD-770-D3   |
| ED30R | DD-701-D3       | DD-770R-D3  |
| ED50  | DD-700-D5       | DD-770-D5   |
| ED50R | DD-701-D5       | DD-770R-D5  |
| ED58  | DD-700-VASC     | --          |
| ED80  | DD-700-D8       | DD-770-D8   |
| ED80R | DD-701-D8       | DD-770R-D8  |



**Pocket-Dop II**

| Natus | DigiDop Classic | DigiDop II |
|-------|-----------------|------------|
| P200  | DD-301-D2       | DD-330R-D2 |
| P300  | DD-301-D3       | DD-330R-D3 |
| P500  | DD-301-D5       | DD-330R-D5 |
| P800  | DD-301-D8       | DD-330R-D8 |



**IMEXDOP CT+**

| Natus | DigiDop Classic | DigiDop II |
|-------|-----------------|------------|
| CT20+ | DD-901-D2       | --         |
| CT30+ | DD-901-D3       | --         |
| CT50+ | DD-901-D5       | --         |
| CT80+ | DD-901-D8       | --         |

