


Newman Medical's Doppler History Timeline



DigiDop

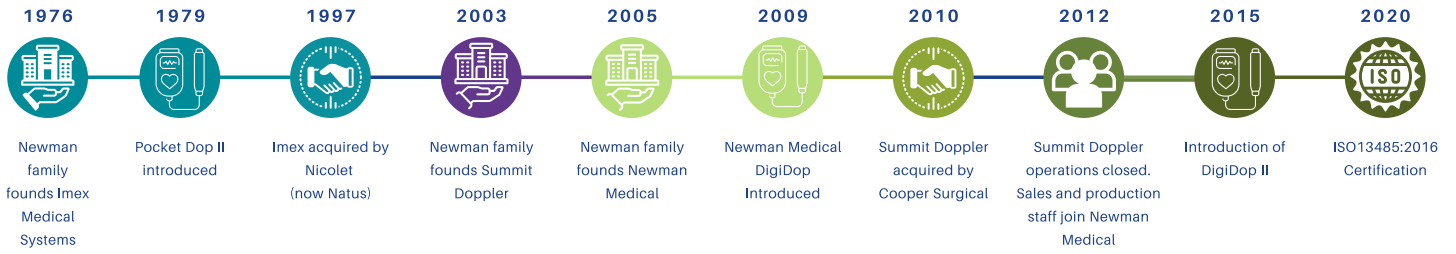

Includes labor, parts, and accidental breakage



American-Made Quality



Includes cords and probes



Waterproof Option

NON-DISPLAY
HANDHELD DOPPLERS
DISPLAY

CooperSurgical Model	SummitDoppler Model	Newman Medical DigiDop Classic	Newman Medical DigiDop II	
 LifeDop 150/150A/150R	LifeDop 150 DD-300	 DD-300	or	 DD-330
	LifeDop 150A DigiDop Classic version not available	DigiDop Classic version not available	or	 DD-330A
	LifeDop 150R DD-301	 DD-301	or	 DD-330R
 LifeDop 250/250AR/250R	LifeDop 250 DD-700	 DD-700	or	 DD-770
	LifeDop 250AR DigiDop Classic version not available	DigiDop Classic version not available	or	 DD-770AR
	LifeDop 250R DD-701	 DD-701	or	 DD-770R
 LifeDop 350 Tabletop	LifeDop 350 DD-901	 DD-901	or	 DD-901

Rechargeable

Audio recording

PROBE ORDERING INFORMATION ON BACK SIDE



Newman Medical Replacement Probe Cross Reference

Vascular Probes



Model D5
5MHz Vascular Probe



Model D8
8MHz Vascular Probe



Model DPPG
Audio TBI PPG Probe

Fetal Probes



Model D2
2MHz Fetal Probe



Model D2W
2MHz Waterproof Fetal Probe



Model D3
3MHz Fetal Probe



Model D3W
3MHz Waterproof Fetal Probe

Professional and Sensitive Dopplers

Quality

- Better sound - Digital Sound Clarity™ produces up to 6x signal of other Dopplers
- Better warranty - 5 year honest warranty all parts & labor. Includes rechargeable batteries on rechargeable units.
- Better support & service - knowledgeable staff and guaranteed service

Sustainability

- Smart Recharge System™ (SRS)
- Significant cost savings
- New rechargeable batteries hold charge for years
- May be recharged more than 2000 times

Innovation

- Continual product development and research

DigiDop Classic



- Standard size
- Longer battery life than similar Cooper products

DigiDop II



- Smaller size
- Ergonomic - fits comfortably in your hand
- Similar battery life to Cooper

DigiDop 901



- Ideal for hospital use
- Great room-filling sound from two speakers on top
- Two probe cradles
- Standard drop-on charging pyramid
- Optional wall mount recharger

How to build your Doppler

DOPPLER CODE

PROBE CODE

Doppler

DD-300	Basic
DD-301	Basic + Recharge
DD-700	HR Display
DD-701	HR Display + Recharge
DD-901	Tabletop + HR Display + Recharge

DigiDop "Classic"

DD-330	Basic
DD-330R	Basic + Recharge
DD-770	HR Display
DD-770R	HR Display + Recharge

DigiDop II

Probe

D2	2MHz fetal probe
D2W	2MHz waterproof fetal probe
D3	3MHz fetal probe
D3W	3MHz waterproof fetal probe
D5	5MHz vascular probe
D8	8MHz vascular probe
DPPG	Audio PPG probe - includes toe cuff

Interchangeable probes can be used on any DigiDop

Vascular



- Great sensitivity
- Ergonomic shape easy to hold, place in sterile glove finger for use in sterile fields
- D8 - 8 MHz** - best for shallow vessels like radial, ankle, foot, digits, penile
- D5 - 5 MHz** - wider beam, deeper penetration. Good for brachial, carotid, femoral, deep veins

Fetal



- Great sensitivity
- Smooth shape, comfortable for patient, easy to hold and clean
- D3 - 3 MHz** - wide angle, ideal for early fetal hearts
- D3W** - a unique 3 MHz waterproof probe
- D2 - 2MHz** - deeper penetration
- D2W** - 2 MHz - waterproof

PPG



- DPPG**
- Unique audio PPG probe for use on digits. Great for toe-brachial studies

Cooper Surgical/Summit Doppler - Newman Medical Part Number Cross-Reference

LifeDop 150/150A/150R

Cooper	DigiDop Classic	DigiDop II
L150-SD2	DD-300-D2	DD-330-D2
L150-SD3	DD-300-D3	DD-330-D3
L150-SD4	DD-300-D5	DD-330-D5
L150-SD5	DD-300-D5	DD-330-D5
L150-SD8	DD-300-D8	DD-330-D8
L150-SDW	DD-300-D2W	DD-330-D2W
L150A-SD2	N/A	DD-330A-D2
L150A-SD3	N/A	DD-330A-D3
L150A-SD4	N/A	DD-330A-D5
L150A-SD5	N/A	DD-330A-D5
L150A-SD8	N/A	DD-330A-D8
L150A-SDW	N/A	DD-330A-D2W
L150R-SD2	DD-301-D2	DD-330R-D2
L150R-SD3	DD-301-D3	DD-330R-D3
L150R-SD4	DD-301-D5	DD-330R-D5
L150R-SD5	DD-301-D5	DD-330R-D5
L150R-SD8	DD-301-D8	DD-330R-D8
L150R-SDW	DD-301-D2W	DD-330R-D2W

*DISCONTINUED SKU

LifeDop 250/250AR/250R

Cooper	DigiDop Classic	DigiDop II
L250-SD2	DD-700-D2	DD-770-D2
L250-SD3	DD-700-D3	DD-770-D3
L250-SD4	DD-700-D5	DD-770-D5
L250-SD5	DD-700-D5	DD-770-D5
L250-SD8	DD-700-D8	DD-770-D8
L250-SDW	DD-700-D2W	DD-770-D2W
L250AR-SD2	N/A	DD-770AR-D2
L250AR-SD3	N/A	DD-770AR-D3
L250AR-SD4	N/A	DD-770AR-D5
L250AR-SD5	N/A	DD-770AR-D5
L250AR-SD8	N/A	DD-770AR-D8
L250AR-SDW	N/A	DD-770AR-D2W
L250R-SD2	DD-701-D2	DD-770R-D2
L250R-SD3	DD-701-D3	DD-770R-D3
L250R-SD4	DD-701-D5	DD-770R-D5
L250R-SD5	DD-701-D5	DD-770R-D5
L250R-SD8	DD-701-D8	DD-770R-D8
L250R-SDW	DD-701-D2W	DD-770R-D2W

*DISCONTINUED SKU

LifeDop 350 Tabletop

Cooper	DigiDop Classic
L350R-SD2	DD-901-D2
L350R-SD3	DD-901-D3
L350R-SD4	DD-901-D5
L350R-SD5	DD-901-D5
L350R-SD8	DD-901-D8
L350R-SDW	DD-901-D2W
*DISCONTINUED SKU	

MADE IN
USA



Scan to learn more about
DigiDop Dopplers



*Crossed out products were discontinued by Cooper Surgical in a letter from Cooper Surgical dated 11/24/23 (ref: C-GL-CSI-000003)