

# DIGIDOP II WITH AUDIO RECORDING Newman Medical



American-Made  
Quality




Includes labor, parts,  
and accidental  
breakage



Audio Recorder



Includes cords  
and probes

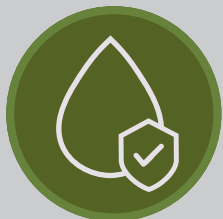

Waterproof fetal  
probe option



## DigiDop

### DigiDop II Audio Recording Models


- DD-770AR
- DD-770A
- DD-330AR
- DD-330A

#### Included




AUDIO RECORDER

#### Standard Options

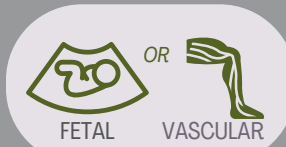


FETAL

VASCULAR

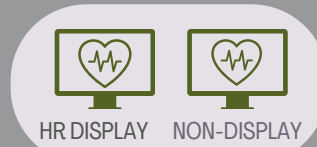

HR DISPLAY

NON-DISPLAY

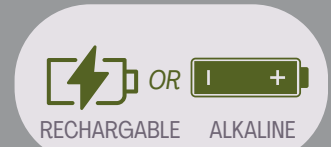

RECHARGABLE

ALKALINE

### Confidently assess blood flow with the DigiDop II Audio Doppler.

Engineered for exceptional audio clarity and ease of use, the DigiDop II empowers healthcare professionals with reliable and accurate vascular and fetal assessments. Featuring advanced technology, a durable design, and an audio recording feature, the DigiDop II is the trusted choice for efficient and effective blood flow analysis in a variety of clinical settings.

### Hear the Difference

- **Digital Sound Clarity (DSC™):** Amplifies the signal up to 6x for exceptional accuracy.
- **Audio Recording:** Capture 30-second audio recordings of vascular or fetal sounds for documentation and analysis.
- **American-made:** Guaranteed quality and craftsmanship.
- **Durable and reliable:** Features an unbreakable probe holder.
- **Interchangeable probes:** For maximum efficiency across various applications.

### Confidence-Inspiring Warranty

- 5-year coverage for all parts and labor, with no hidden fees.
- Rechargeable convenience (on applicable models).

### Newman Medical: Your Ultrasound Partner

- Pioneers in Doppler technology since the 1970s.
- Founders of Imex Medical Systems and Summit Doppler.
- Committed to continuous development of cutting-edge solutions.
- Focused on building long-term partnerships and exceeding customer expectations.

# DigiDop Professional and Sensitive Dopplers



Over **45 years in the marketplace**; a testament to our longevity and success



**Made in the USA** for superior quality and performance



**5-year all-inclusive warranty** for maximum protection and minimized downtime, headaches



**Best-in-class customer support**; we always answer the phone

## DigiDop Classic

- Standard size
- Superior battery endurance compared to competitors



## DigiDop II

- Compact size
- Comfortable palm-fitting design
- Equivalent battery endurance to competitors
- Audio recorder option



## DigiDop 901

- Ideal for hospital use
- Clear stereo sound from dual top speakers
- Two probe cradles
- Standard drop-on charging pyramid
- Optional wall mount recharger



## Vascular Probes



**Model D8**  
8MHz Vascular Probe

OR  
**Model D5**  
5MHz Vascular Probe



**Model D8-SD**  
8 MHz SteadyDop Vascular Probe

OR  
**Model D5-SD**  
5 MHz SteadyDop Vascular Probe



**Model DPPG**  
Audio TBI PPG Probe

## Fetal Probes



**Model D2**  
2MHz Fetal Probe

OR

**Model D3**  
3MHz Fetal Probe



**Model D2W**  
2MHz Waterproof Fetal Probe

OR

**Model D3W**  
3MHz Waterproof Fetal Probe

Interchangeable probes can be used on any DigiDop

## Vascular

- Great sensitivity
- Ergonomic shape easy to hold, place in sterile glove finger for use in sterile fields
  - **D8 - 8 MHz** - best for shallow vessels like radial, ankle, foot, digits, penile
  - **D5 - 5 MHz** - wider beam, deeper penetration. Good for brachial, carotid, femoral, deep veins

## Fetal

- Great sensitivity
- Smooth shape, comfortable for patient, easy to hold and clean
  - **D3 - 3 MHz** - wide angle, ideal for early fetal hearts
  - **D3W** - a unique 3 MHz waterproof probe
  - **D2 - 2 MHz** - deeper penetration
  - **D2W** - 2 MHz - waterproof

## DPPG

- Unique audio PPG probe for use on digits. Great for toe-brachial studies

## How to build your Doppler

DOPPLER CODE

PROBE CODE

### Doppler

DD-300	Basic
DD-301	Basic + Recharge
DD-700	HR Display
DD-701	HR Display + Recharge
DD-901	Tabletop + HR Display + Recharge

DigiDop II

DD-330	Basic
DD-330R	Basic + Recharge
DD-770	HR Display
DD-770R	HR Display + Recharge

### Probe

D2	2MHz fetal probe
D2W	2MHz waterproof fetal probe
D3	3MHz fetal probe
D3W	3MHz waterproof fetal probe
D5	5MHz vascular probe
D8	8MHz vascular probe
DPPG	Audio PPG probe - includes toe cuff